



## ACRO™ RSV/Influenza A+B Combo Rapid Test Cassette

A rapid test for the qualitative detection of Respiratory Syncytial Virus Antigen and influenza A and B antigens in Nasopharyngeal swab or Nasal aspirate specimens.

### Reliable Rapid Tests for

- Respiratory Syncytial Virus (RSV)
- Influenza A+B


### Features & Benefits

- Results in 15 minutes
- Easy to use
- Sample types:
  - ▶ Nasopharyngeal swab
  - ▶ Throat swab
  - ▶ Nasal aspirate specimen

### Benefits of Rapid Diagnosis

- Guide appropriate therapy
- Avoid unnecessary antibiotic use
- Reduce cost of treatment


### Direction for use

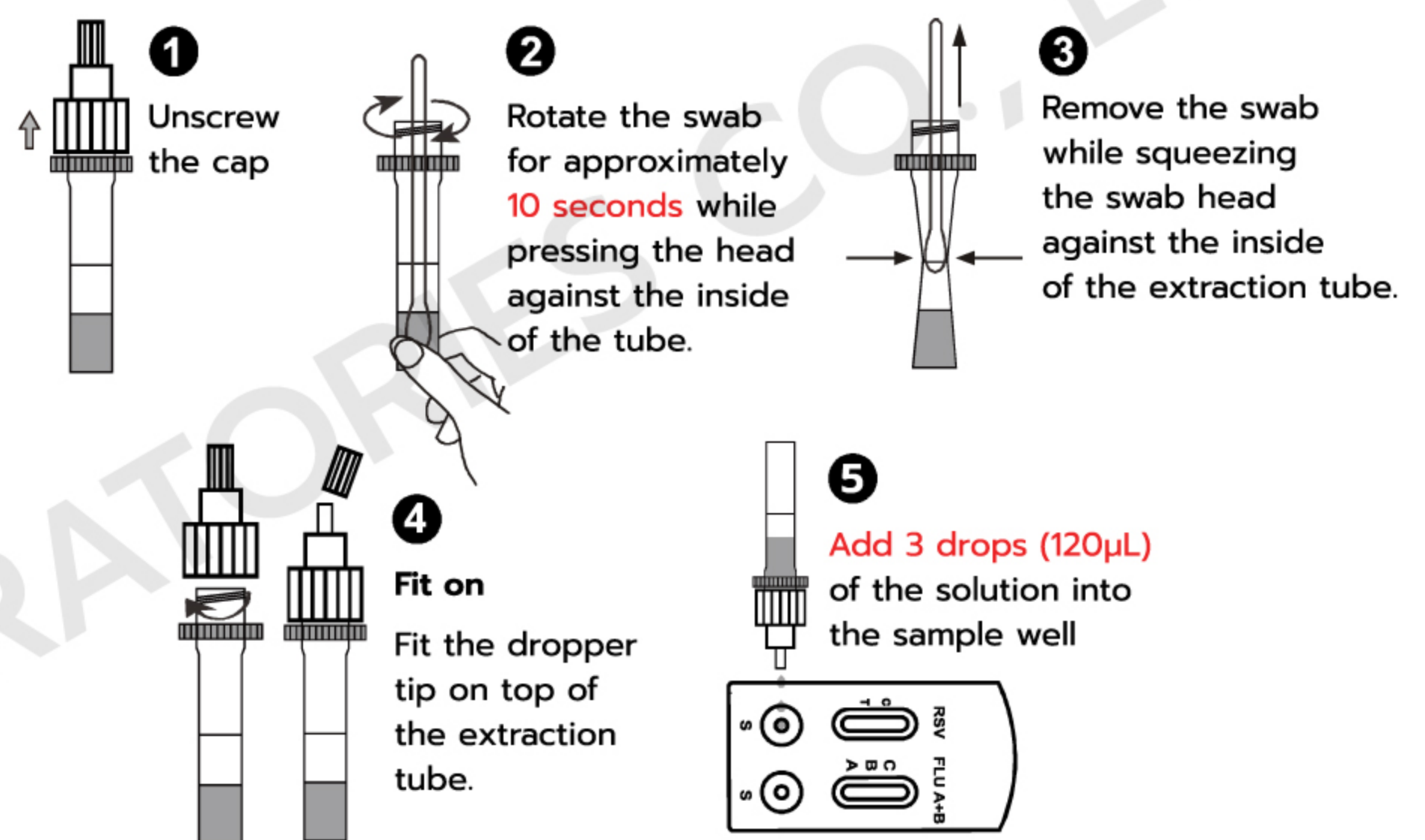

Read the result at 15 minutes

**Note: Do not read result after 20 minute**

### Interpretation of Results

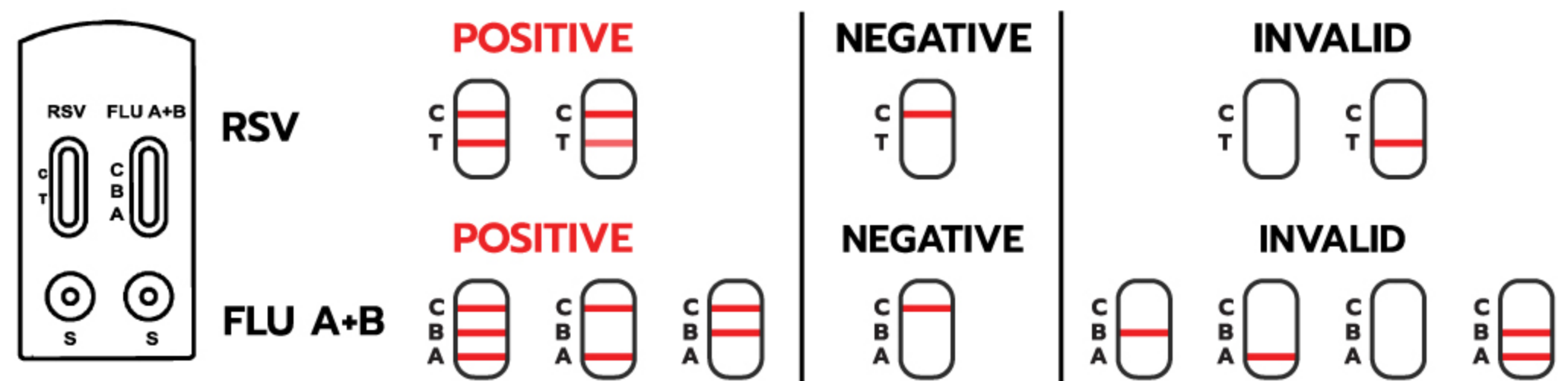

**POSITIVE:** RSV Two colored lines appear  
Influenza A or B-Two distinct colored lines appear.  
Influenza A and B-Three distinct colored lines appear.

**NEGATIVE:** One colored line appears in the control line region (C).

**INVALID:** No line appears in the control line region (C).

ACRO™ provides the rapid tests with high accuracy

Analyte	Sensitivity	Specificity	Accuracy
Flu A	99.0%	98.9%	98.9%
Flu B	97.7%	99.0%	98.6%
RSA	92.7%	98.0%	95.6%

See more about ACRO™ RSV/Influenza A&B performance in Package insert.